




SOMNODENT PRODUCT CHART

	SIGNATURE LINE			STANDARD LINE		
						
	FUSION	FLEX	CLASSIC	AIR/AIR+	HERBST ADVANCE	MORNING REPOSITIONER
WARRANTY	3 YEARS	3 YEARS	3 YEARS	1 YEAR	2 YEARS	1 YEAR
PRODUCT FEATURES	Rapid advancement with interchangeable wings providing an 8.5 mm range and fine screw adjustments at 0.1 mm increments	Features a soft inner liner made with SMH BFlex material, a proprietary soft polymer that is molecularly bonded to the acrylic device	The original innovation in customized oral devices, the Classic is held in place using medical grade stainless steel ball clasp	Available with ball clasps or standard Poly+ liner for retention. Affordable price, optimizes low reimbursement scenarios	8 mm range of calibration with a visual calibration indicator	Hands-free bite deprogrammer incorporating raised blocks designed to appropriately reduce or eliminate the effects of bite change
RECOMMENDED CLEANING AND STORAGE	Rinse with cold or lukewarm water, soak with a SomTab in a container of water, rinse, then store in a container of clean water	Rinse with cold or lukewarm water, soak with a SomTab in a container of water, rinse, then store in a container of clean water	Rinse with cold or lukewarm water, soak with a SomTab in a container of water, rinse, then store in a container of clean water	Rinse with cold or lukewarm water, soak with a SomTab in a container of water, rinse, with lukewarm water, dry and store in a dry container	Rinse with cold or lukewarm water, soak with a SomTab in a container of water, rinse, then store in a container of clean water	Rinse with cold or lukewarm water, dry and store dry
TURNAROUND TIME	14 WORKING DAYS	14 WORKING DAYS	14 WORKING DAYS	14 WORKING DAYS	14 WORKING DAYS	14 WORKING DAYS

PATIENT CONDITION	SUGGESTED FUNCTION	SOMNOMED ORAL DEVICE RECOMMENDATIONS					
		✔ Device is specially made for this patient type	✔ Device will work very well for this patient type	✔ Device will work for this patient type	✔ Device will work for this patient type	✔ Device will work for this patient type	✔ Device will work for this patient type
Open and closing mouth is a priority for the patient	Allows mouth to open/close more freely without metal hinges	✔	✔				✔
Larger Tongue	Increased tongue space for patient					✔	✔
Bruxism/TMD Sufferer	Increased tongue space for patient					✔	✔
Smaller Mouth	Shorter adjustment arms for compatability / Fusion is slimmer with less bulk due to prefabricated wings	✔				✔	✔
Smaller Mouth & Bruxism/TMD Sufferer	Shorter adjustment arms for compatability and durability for clenching & grinding					✔	✔
Mouth Breathers	Option of anterior opening for capability of mouth breathing while device is worn	✔				✔	✔
Extreme vertical opening / Patient with large vertical	Option of elastic retention (ER) hooks	✔	✔			✔	✔
Transportation patient / Positional Therapy Patient	DentiTrac not only allows for compliance tracking, but sleep position as well	✔	✔			✔	✔

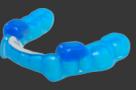